



Harvard ZirconCore

**Dual cure composite with nano zirconium dioxide.
For core build-up and post cementation.**


- Excellent consistencies for both indications:
post cementation and **core build-up**
 - Optimal thixotropic properties
 - Optimal hardness: **“Cuts like dentin”**.
Avoids grooves from circular grinding
 - Very good mechanical properties for permanent restorations
 - Contains a substantial amount of Nano-zirconia
 - Dual-curing
 - Fluoride release
 - Very good radiopacity
 - Self-mixing material with the Minimix syringe
- For optimal results use **Harvard ZirconCore** with the adhesive **Harvard Bond SE Dual**

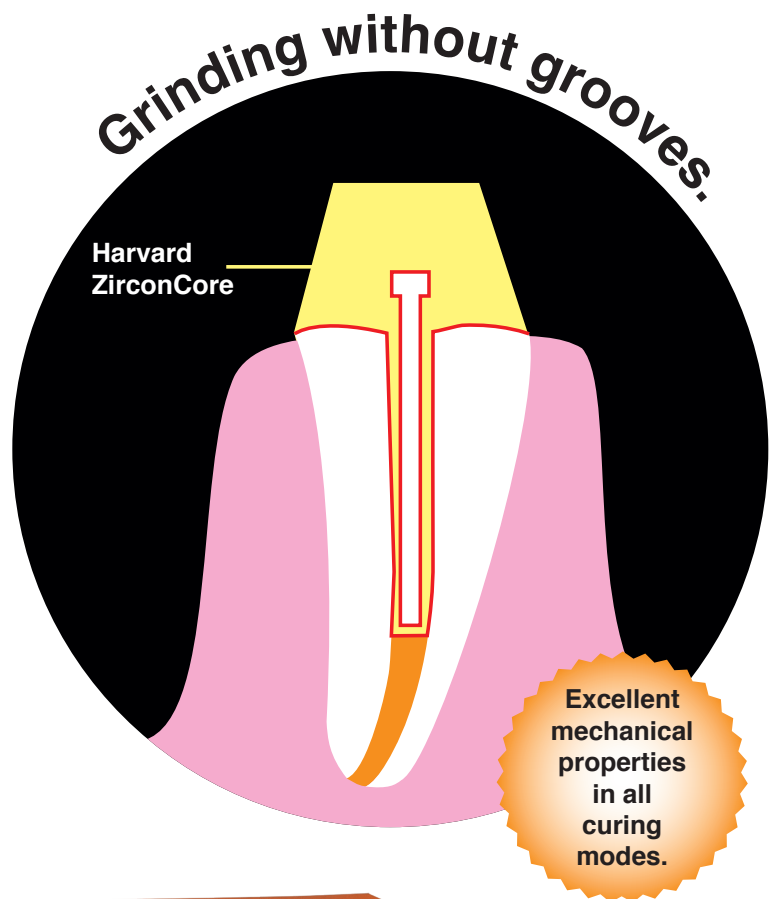

Physical Properties

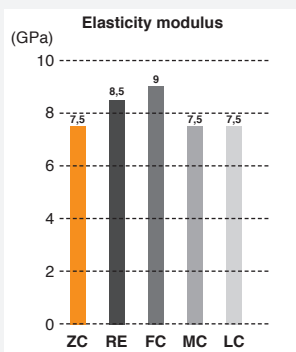


(light cure resp. dual cure)

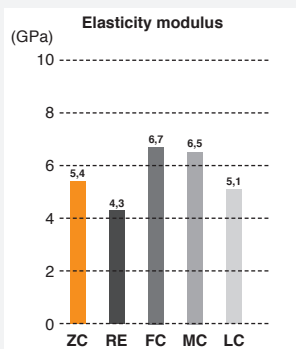

(dark cure)

ZC Harvard ZirconCore
 RE Voco Rebuilda
 FC Dentsply FluoroCore 2+
 MC Ivoclar Vivadent MultiCore Flow
 LC DMG LuxaCore Z Dual Automix

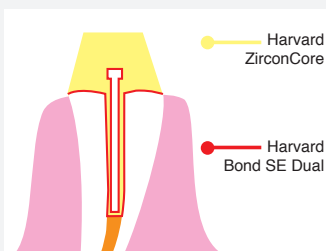

Post cementation and core build-up



1. Prepared and bonded root canal
2. Application into the root canal
3. Placement of the post
4. Light cured core build-up
5. Shaped core build-up

Technical data

Working time (23 °C / 73 °F)	1:30 min
Setting time (self cure)	3:30 min
Light cure (optional)	40 sec

Article information

Article information	Order no.
Harvard ZirconCore 5 ml minimix syringe, 10 mixing tips, 10 intra tips long, Shade A3	7083600
Harvard Mini 1:1 O-Brown refill bag with 50 mixing tips	7083610
Harvard IntraTips long refill bag with 50 intra tips long	7083620
Harvard Bond SE Dual 2 x 5 ml bottle, 50 Microbrush®, 5 mixing wells	7083602
Harvard Bond SE Dual Refill 2 x 5 ml bottle	7083605

Liability is excluded for all printing errors and omissions.
 Before using our Harvard products, the respective directions for use should be observed in every case.
 Internal data by Harvard Dental International.

Harvard Distribution Partner.

We recommend
Harvard Bond SE Dual:
 the dual cure,
 self etch
 adhesive



125
 JAHRE
 1892-2017
 HARVARD



HARVARD®

Marke und Qualität seit 1892

Harvard Dental International GmbH
 Margaretenstr. 2 - 4, 15366 Hoppegarten, Germany
 Phone: + 49 (0) 30/99 28 978-0
 Fax: + 49 (0) 30/99 28 978-19
 info@harvard-dental-international.de